


Installing HF108B Filter Wheel onto Evident Microscopes

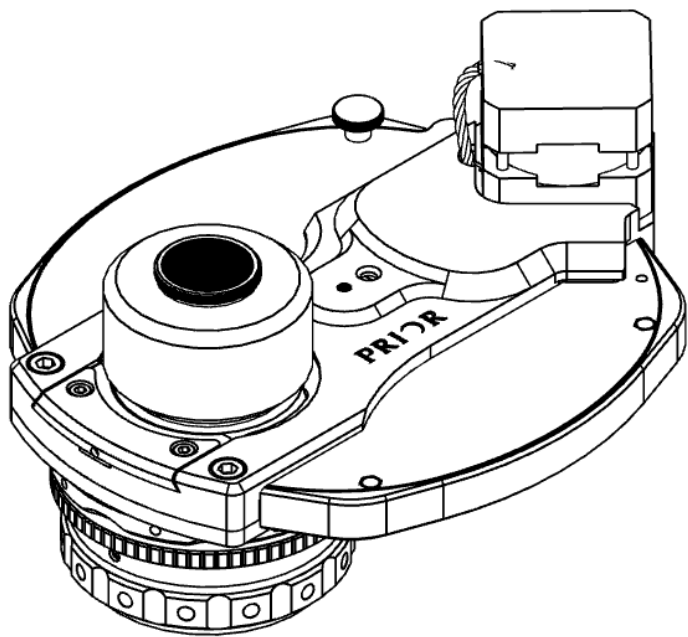


Contents

Installation onto Evident Microscope	2
Changing the installed camera mount	7
Making parfocal adjustments.....	11
Changing camera orientation.....	14
Installing and removing filters	16
About Prior Scientific.....	19
Our mission	19
Technical Support.....	19


Installation onto Evident Inverted Microscope



To install the filter wheel on the left side port, HF108B models with the suffix EV will be required, and for the right side port, the HF108B models with the suffix IX85RH is required.

1. Loosen grub screw on the microscope side port.

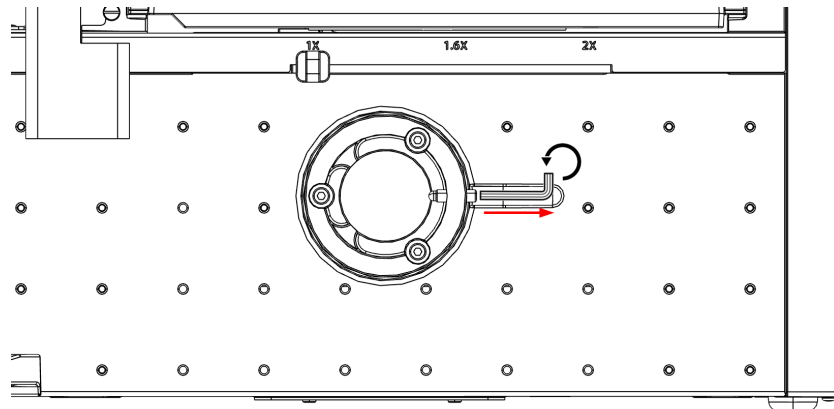


2. Place filter wheel dovetail into microscope receiving dovetail.

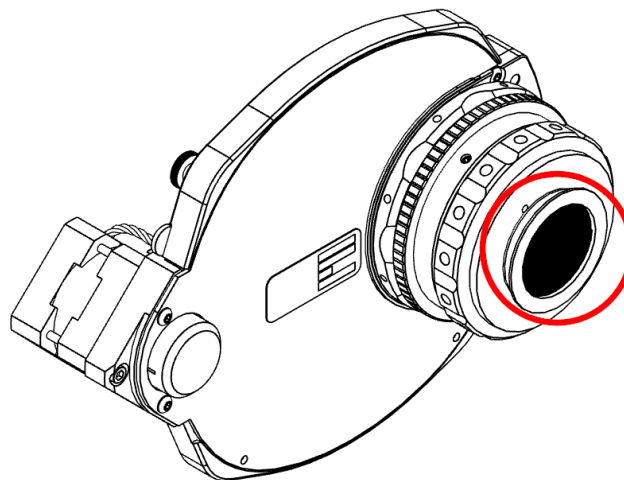

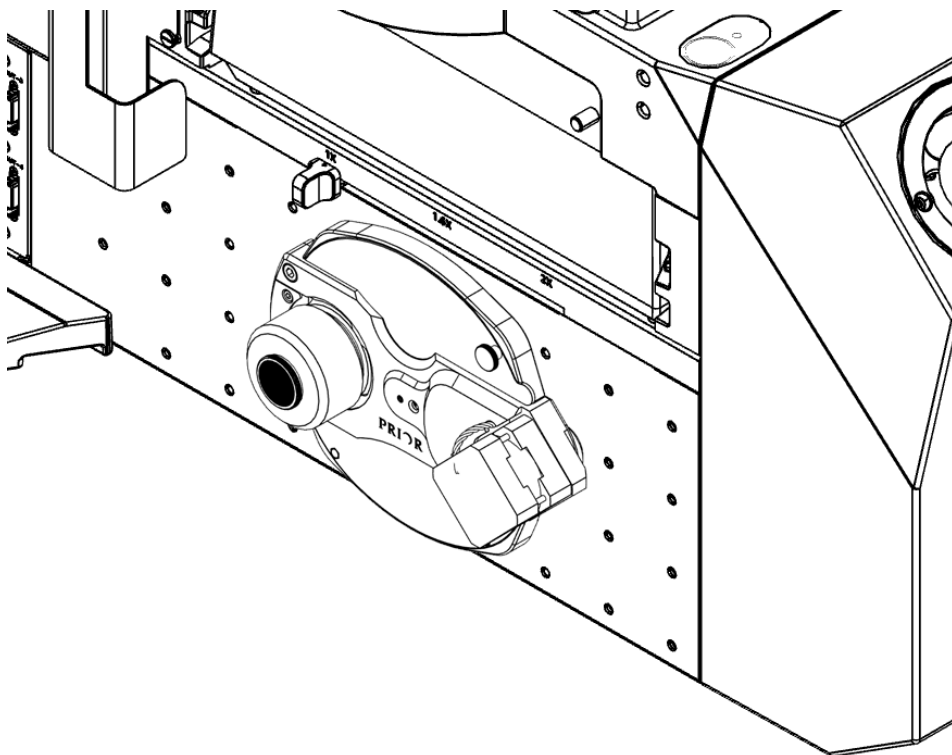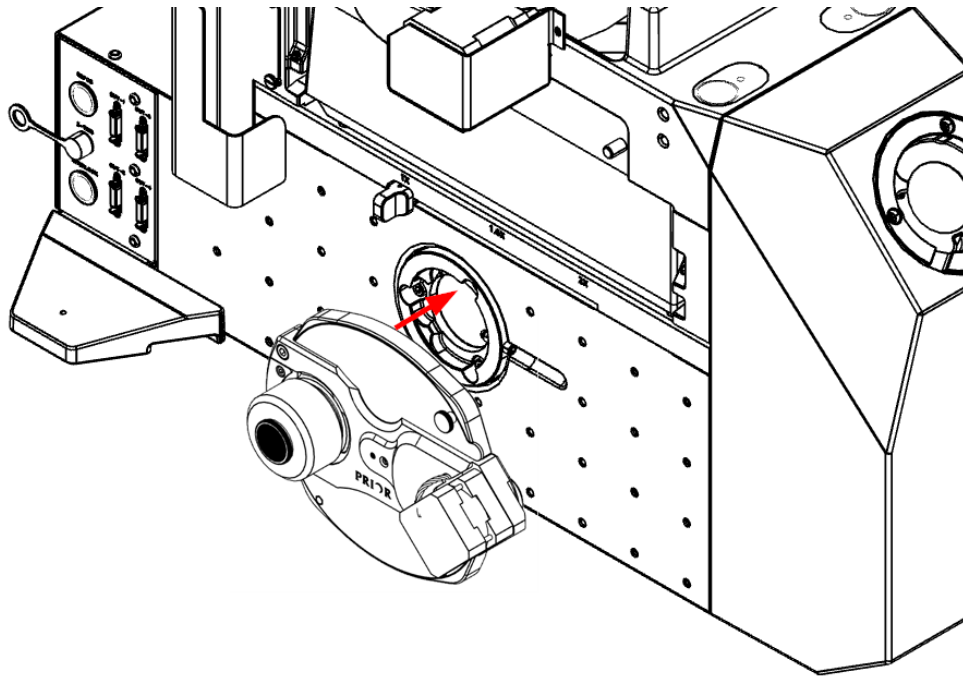


3. Tighten the grub screw to secure filter wheel in place.

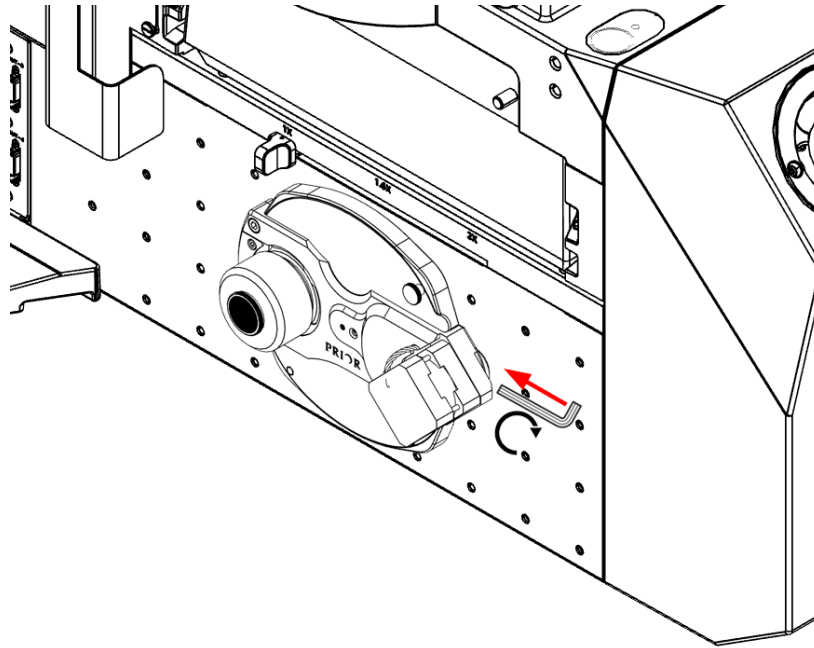



To Install onto Evident Upright Microscope

1. HF108B models with the suffix -EV can also be installed onto an upright microscope that is compatible with the C-mount adapter part combination U-TV1X plus U-CMAD3. To allow fitting, first remove any cameras that may be attached, and loosen the grub screw indicated below.

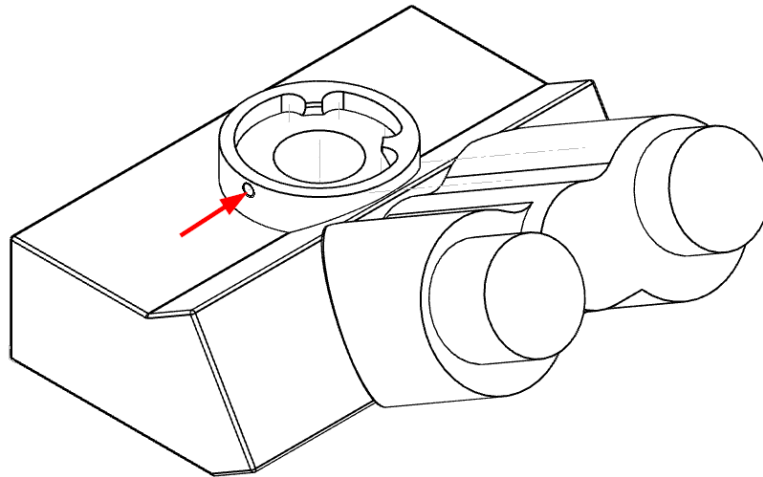


2. Align the microscope adapter side of the filter wheel with the trinocular head.

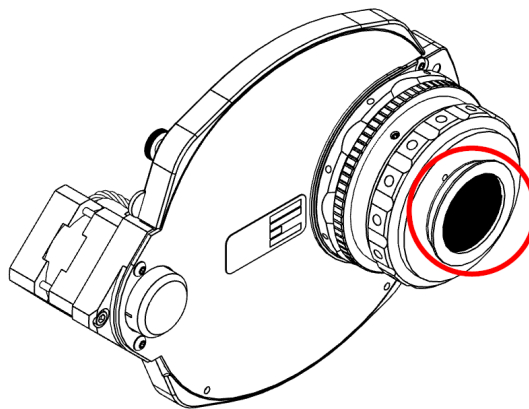

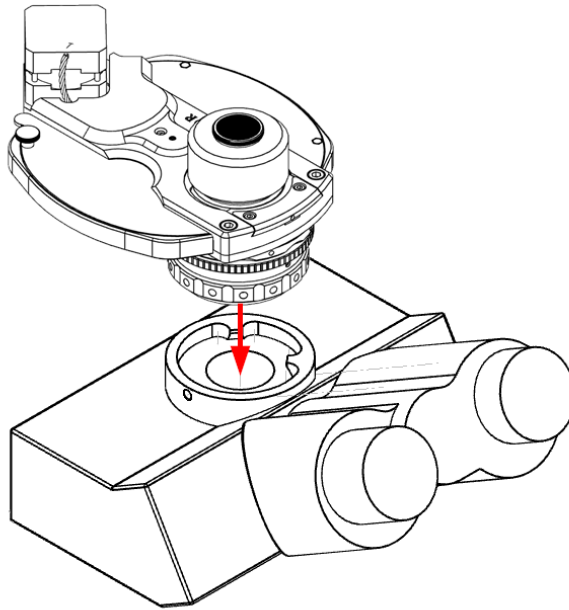


- 3.** Secure the filter wheel into place with the grub screw.

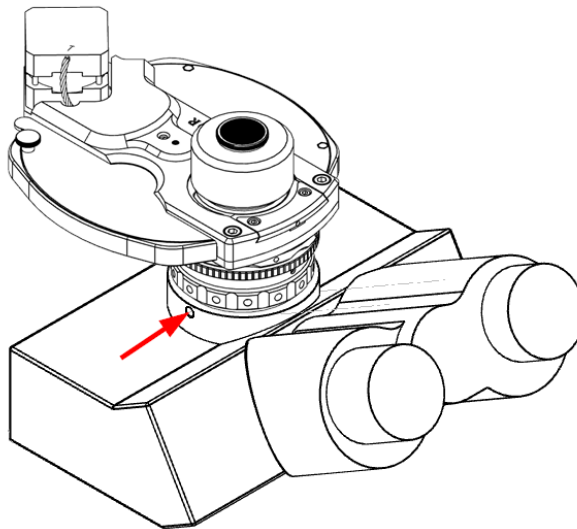


Changing the installed camera mount


1. Unscrew the two bolts securing the camera mount.

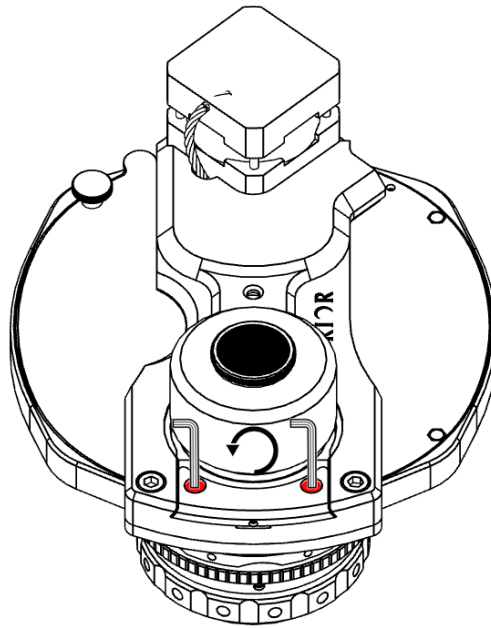


2. Slide out the camera mount locking bar.

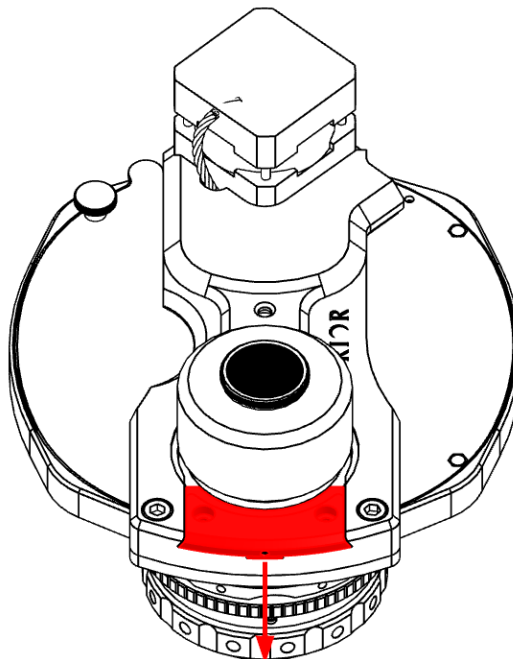


3. Slide out the brass dovetailed piece, noting the correct orientation.

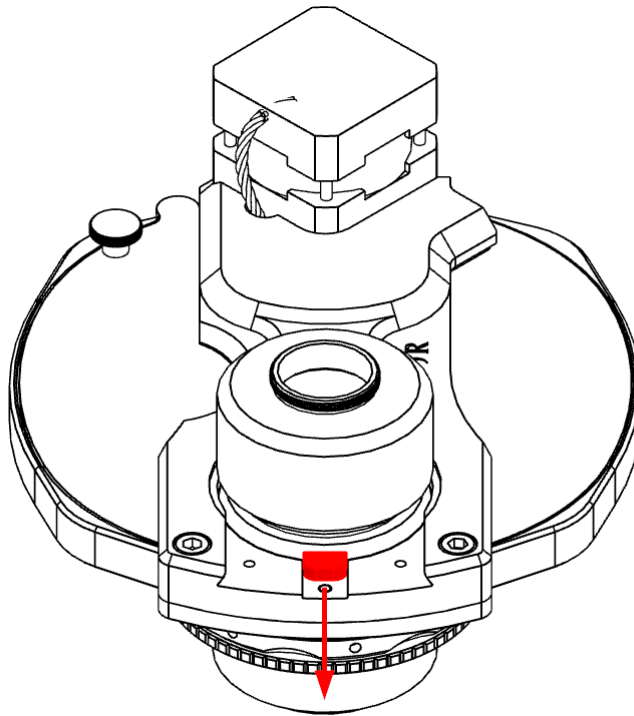


4. Slide out the previous camera mount.

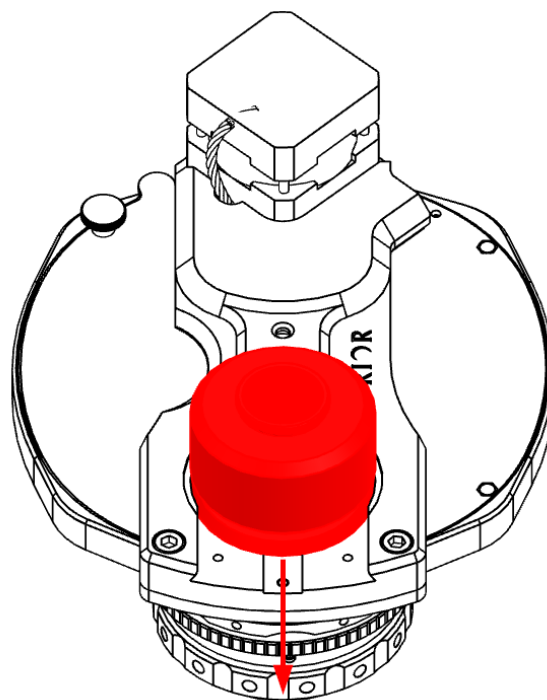


5. Slide in the new camera mount

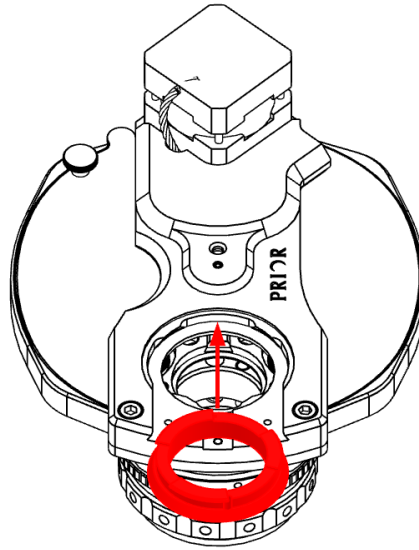


6. Replace the brass dovetailed piece, maintaining the previous orientation.

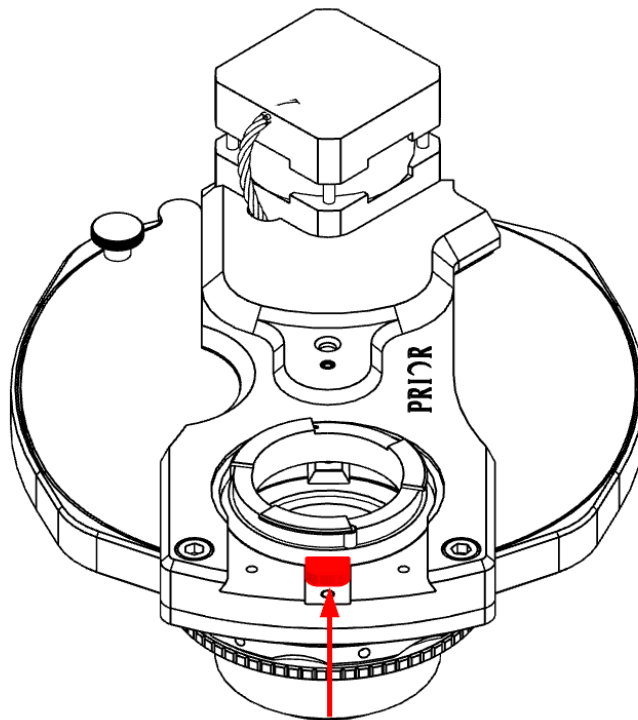

7. Replace the camera mount locking bar.

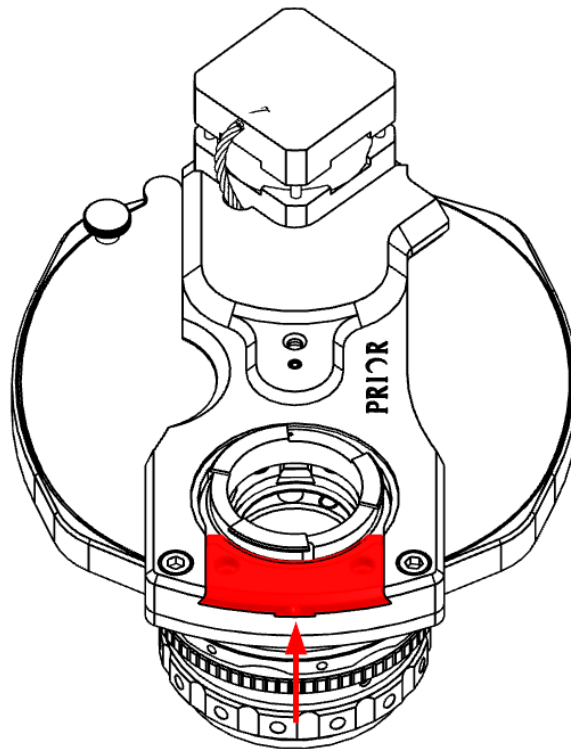



8. Re-screw the locking bar into place.

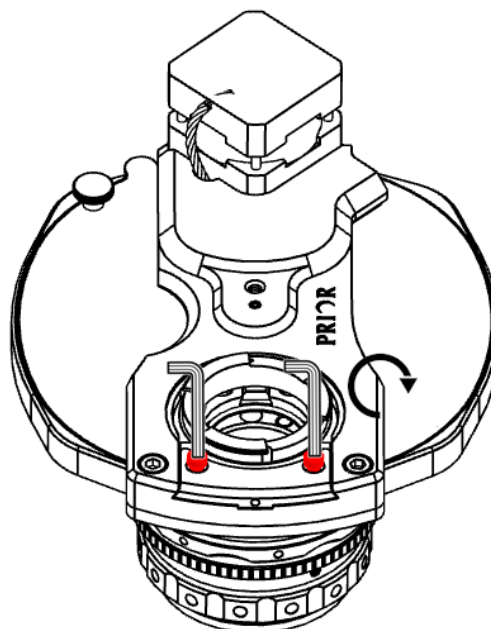


9. Replacement is now complete.


Making parfocal adjustments

1. Unlock the parfocal adjuster by releasing the three tightened screws, indicated in the images below.

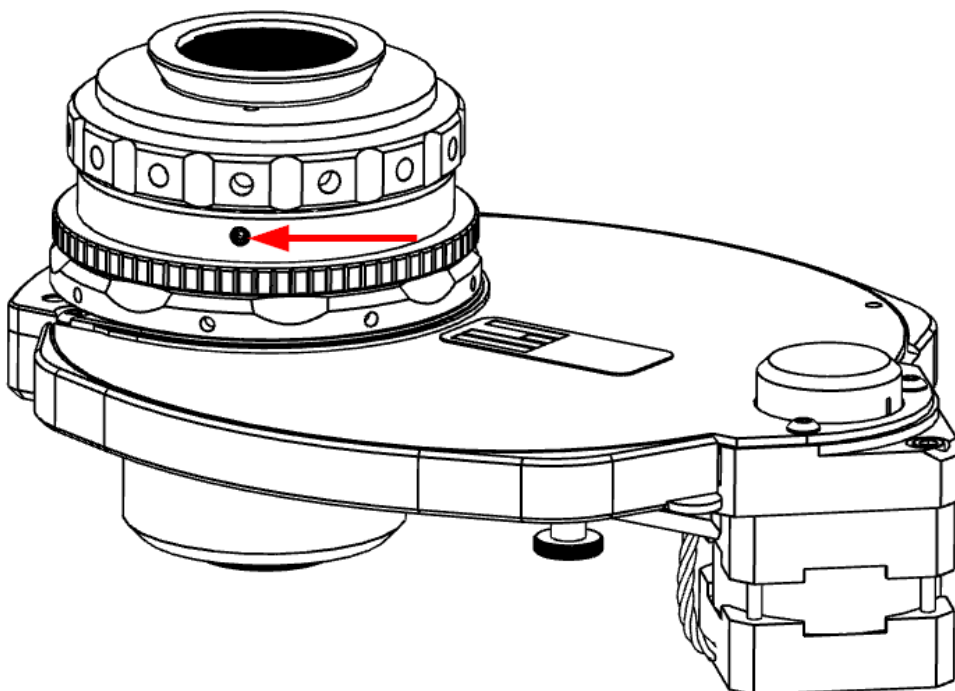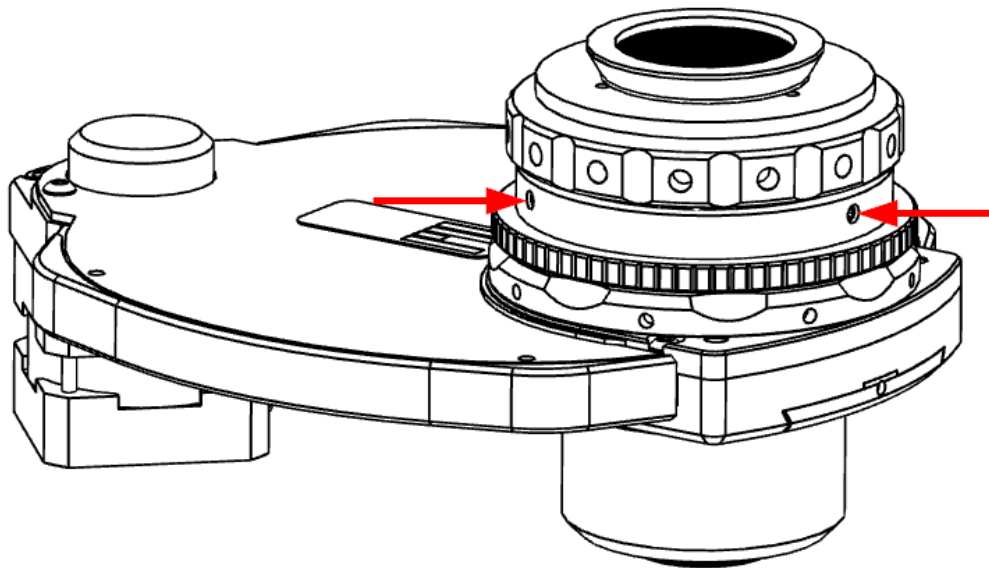


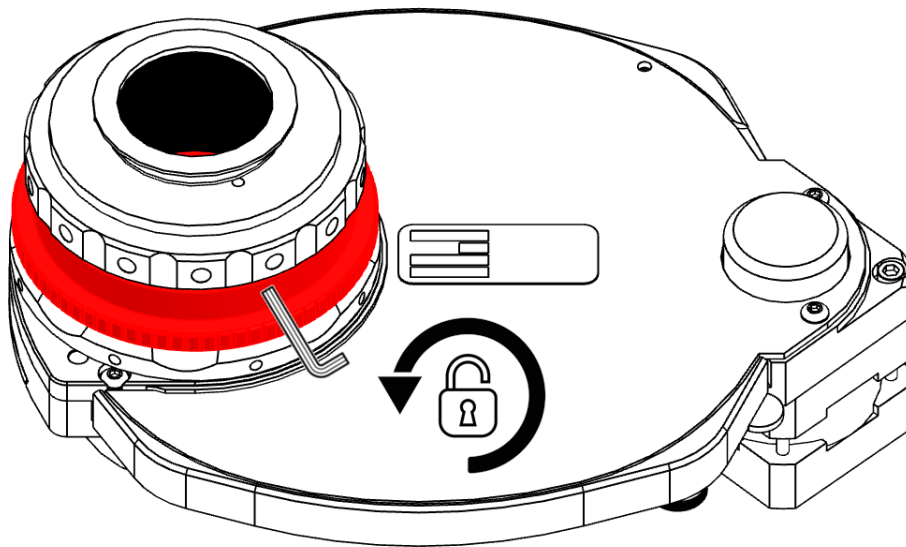

2. Adjust the adapter focus ring until parfocality is achieved. A hex key may be required to reach the focus ring when the filter wheel is mounted in position.
3. Turning the adapter anticlockwise will extend the path length, turning the adapter clockwise will reduce path length.

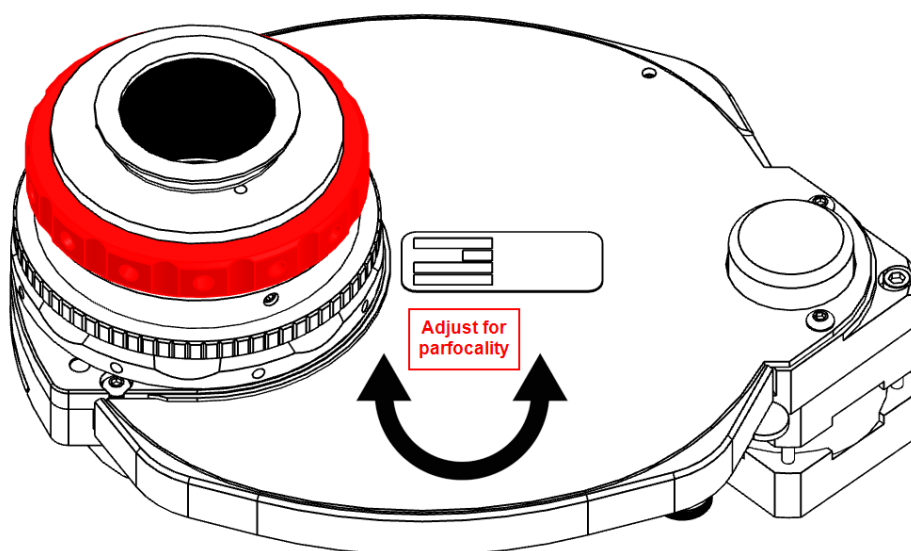


4. Once parfocal, re-tighten all available screws

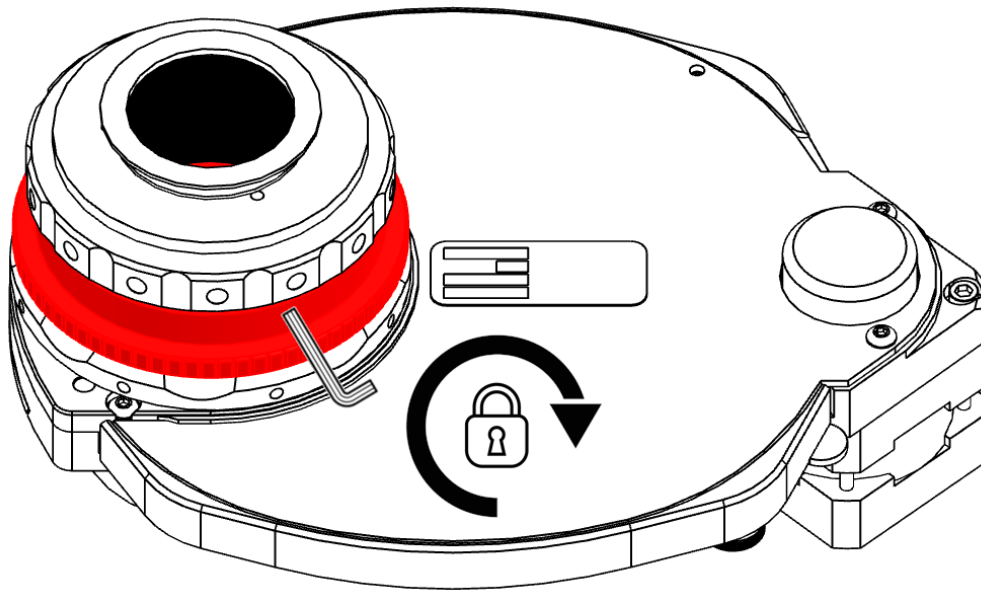


5. Parfocality adjustment is now complete.

Changing camera orientation


1. Loosen circled grub screw.

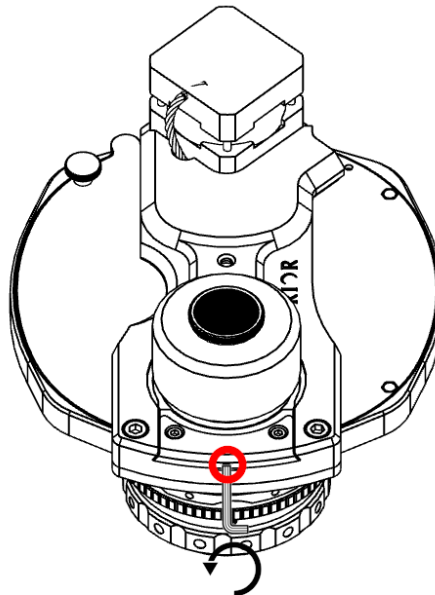

2. Rotate camera mount to desired position.

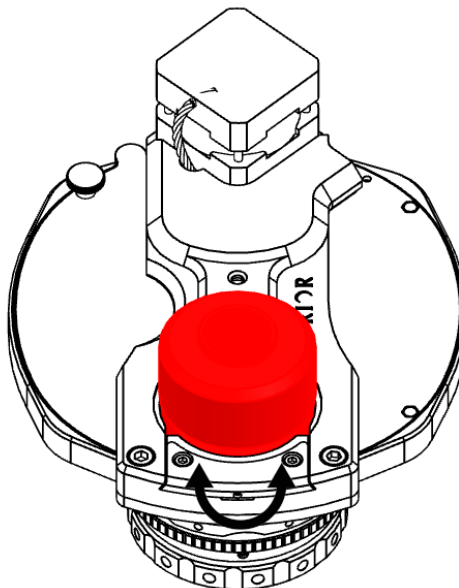



3. Re-tighten circled grub screw.

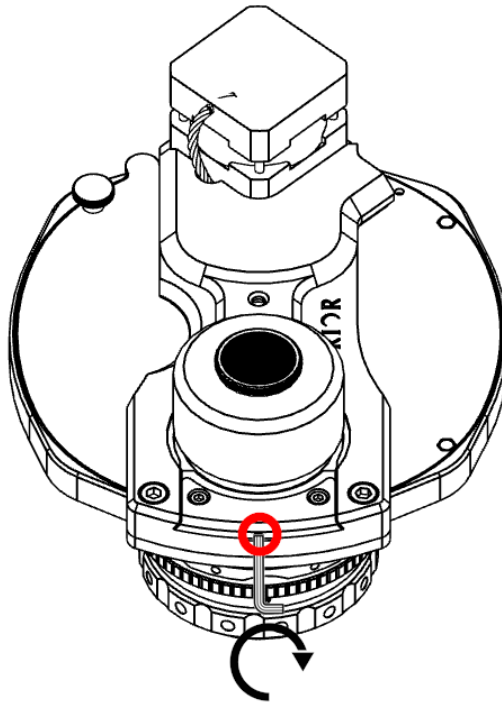


4. Camera orientation change is complete.

Installing and removing filters


1. Remove the filter cover.

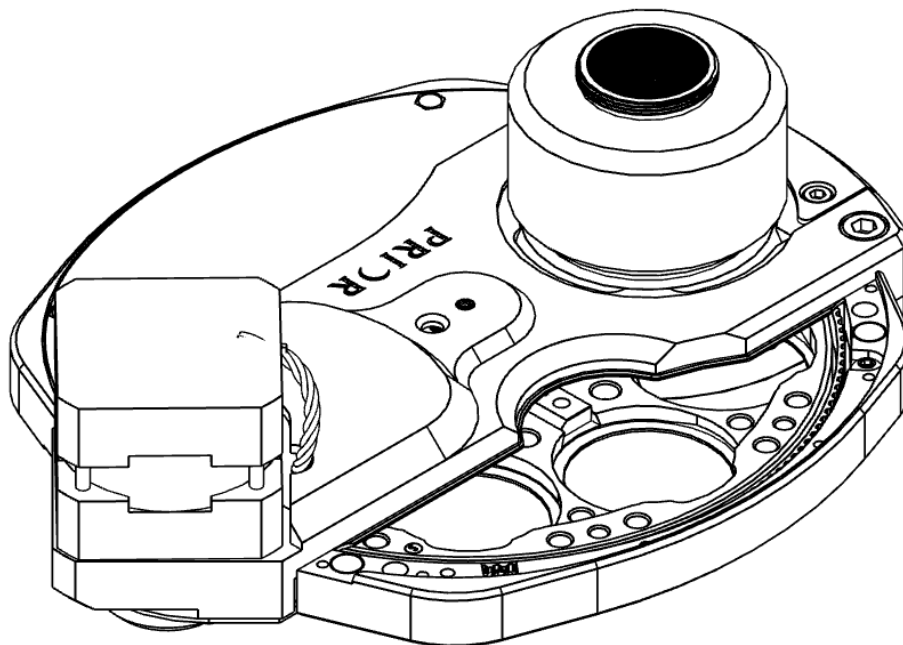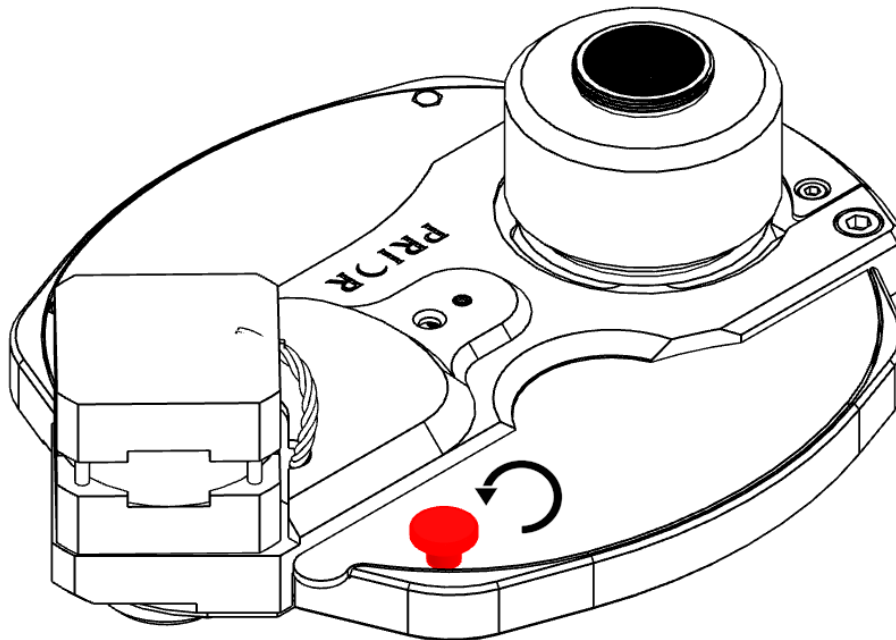


2. Confirm filter position.

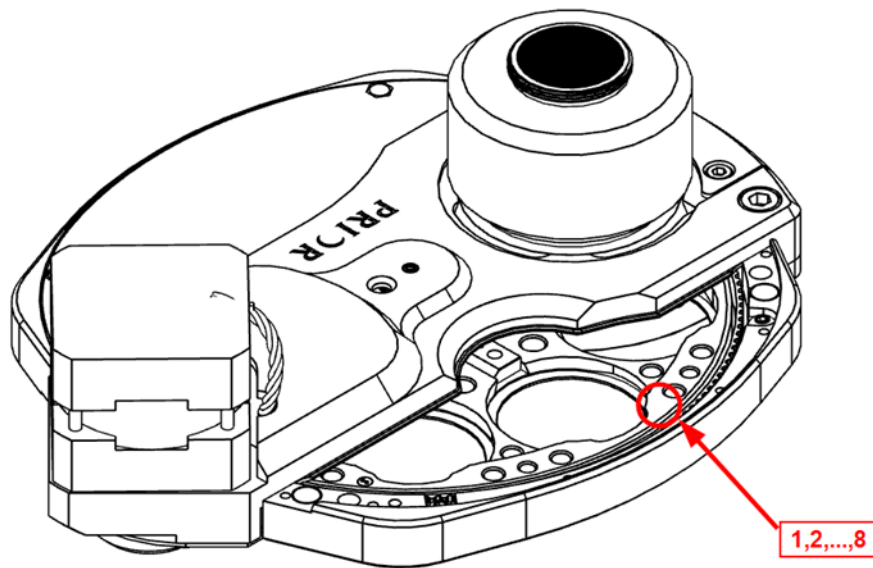


3. Insert filter into wheel by pulling back retainer with hex key.

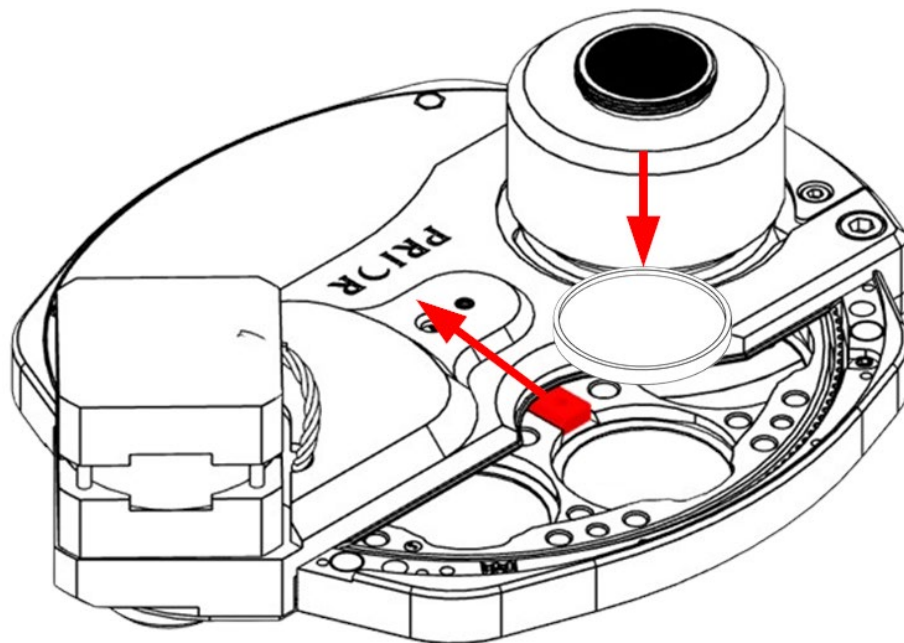

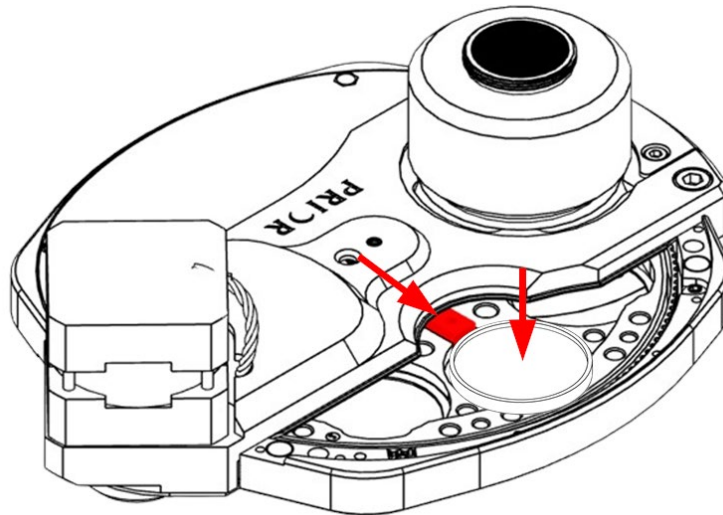

4. Repeat for all desired filters.

5. Reattach filter cover.

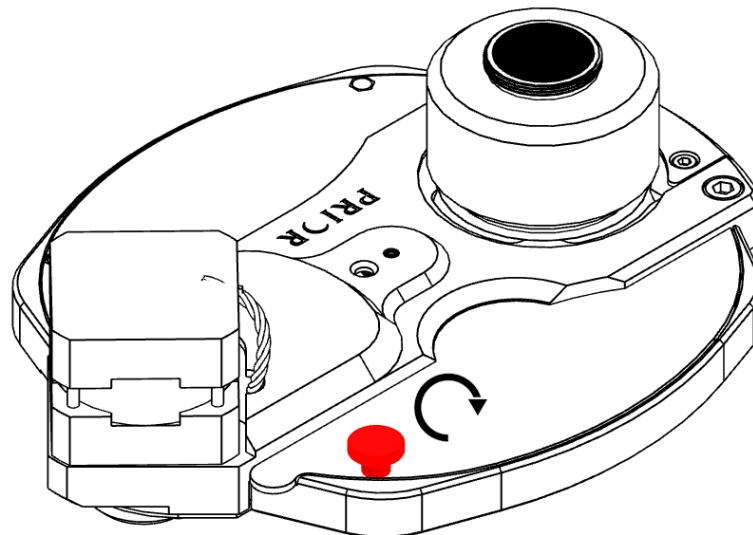

6. Filter installation is now complete. Filter wheel should be rehommed before use. This may be autohomed if setup previously with the 7,1,A serial command or via the serial command 7,1,H (where 1 is the filter port on the Prior Controller).

7. To remove filters, pull back retainer with hex key and lift filter from cell.

About Prior Scientific

Prior Scientific has manufactured precision optical and electromechanical equipment and instrumentation for over 100 years.

As the leader in the development of systems for precision motion control and microscope automation, we work closely with our partners and customers to develop customized optical and automation systems as well as providing a wide range of standard microscopy components.

We partner with system integrators and OEMs to develop prototypes, custom-built products, and complete imaging system solutions according to individual customer requirements for a variety of applications. Our integrated design and manufacturing facilities enable us to seamlessly transition from concepts to production systems. Typical applications for our products are found in photonics, life science research, material science, nanotechnology, medical devices, metrology, and semiconductor manufacturing.

Our mission

We understand that every problem requires a unique approach. Our team works closely with customers and partners to identify their specific needs and design a cost-effective and tailored solution that meets their unique requirements. We continue to push scientific boundaries and constantly build on and improve our expertise in automated precision optical and mechanical systems.

Technical Support

If you have a technical question, you can contact our expert support team:

Tel: +44 (0) 1223 881711

Email: techsupportUK@prior.com

Visit prior.com and use our online chat

UNITED KINGDOM

Prior Scientific Instruments Ltd.
Units 3-4 Fielding Industrial Estate
Wilbraham Road, Fulbourn
Cambridge, CB21 5ET
United Kingdom
Email: inquiries@prior.com
Phone: +44 (0)1223 881711

U.S.A.

Prior Scientific, Inc.
80 Reservoir Park Drive
Rockland, MA. 02370
U.S.A.
Email: info@prior.com
Phone: +1 781 878 8442

GERMANY

Prior Scientific Instruments GmbH
Maria-Pawlowna-Str. 4
D-07743, Jena, Germany
Email: jena@prior.com
Phone: +49 (0)3641 242 010

JAPAN

Kayabacho 3rd Nagaoka Bldg 10F,
2-7-10, Nihonbashi Kayabacho, Chuo-Ku,
Tokyo103-0025, Japan
Email: info-japan@prior.com
Phone: +81 (0)3 5652 8831

CHINA

Prior Scientific Instruments (Suzhou) Ltd.
Room 118, Meilihua Hemu Park
No. 393 Suhong Middle Road, Suzhou Industrial Park
Suzhou, 215000, China
Email: info-china@prior.com
Phone: +86 (0)512 6617 5866

