

H116SPN

Motorised semi-conductor wafer inspection shuttle stage for Nikon wafer handling systems



The Prior Scientific motorised shuttle stage is designed to be used for wafer inspection, using the Nikon NWL200 and NWL860 wafer handling systems. The stage also reduces operator error and fatigue whilst increasing the rate at which wafers are inspected.

The H116SPN automatically loads and unloads wafers. When operating the shuttle system automatically moves to the correct load position and activates the loader switch, notifying the loader that the wafer chuck is ready to receiving the wafer. After transfer to the shuttle is retracted allowing for inspection. Either a manual inspection process using a joystick, or a pre programmed inspection pattern, are possible.

The H116SPN has a large travel range coupled with superb repeatability and resolution. When used in conjunction with Prior's ProScan® system, it ensures accuracy, precision and ultimately a more effective and efficient wafer inspection process.


- Travel range of 250 x 210 mm
- Resolution of 0.04 μm
- Repeatability ± 0.2 μm
- · Automatically loads and unloads wafers
- Compatible with the ProScan III system and Nikon NWL200/NWL860 wafer handling systems



H116SPN

Motorised semi-conductor wafer inspection shuttle stage for Nikon wafer handling systems

Dimensions

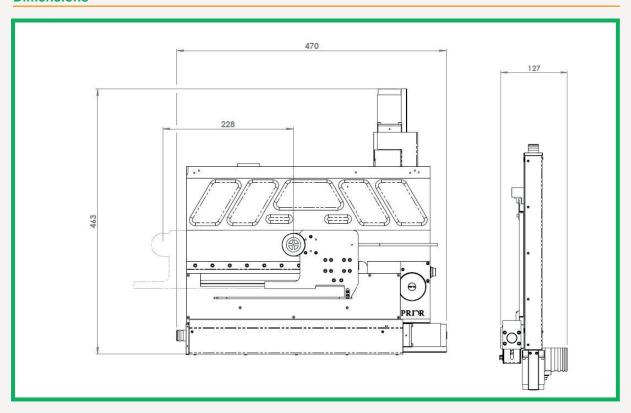

Specification*

| Feature | Specification |
|--------------------------------|---------------|
| Bi-directional repeatability | ± 1.3 μm |
| Uni-directional repeatability | ± 0.2 μm |
| Minimum step size (resolution) | 0.04 μm |
| Recommended speed | 40 mm/s |
| Maximum travel range | 250 x 210 mm |
| Ball Screw pitch | 2 mm |
| Motor Type | 200 SPR |
| Weight | 9 Kg |

^{*} To obtain these values requires the use of a Prior ProScan II or above controller. Using the Prior method of testing with backlash correction enabled. All repeatability is uni-directional.





Worldwide distribution

Prior Scientific Ltd Cambridge, UK T. +44 (0) 1223 881 711 E. uksales@prior.com Prior Scientific Inc Rockland, MA USA T. +1 781-878-8442 E. info@prior.com Prior Scientific GmbH Jena, Germany T. +49 (0) 3641 675 650 E. jena@prior.com Prior Scientific KK Tokyo, Japan T. +81-3-5652-8831 E. info-japan@prior.com

© 2015 Prior Scientific Instruments Ltd. Specifications subject to change at any time. E & O E.